


Tune Check Report

Operator Name admin
Acq/Data Batch D:\Agilent\ICPMH\1\UserTune_8900.b
Acq. Date-Time 24-May-24 12:55:15 PM
Report Comment ---
Instrument Name G3665A SG18461073

[No Gas]

Sensitivity

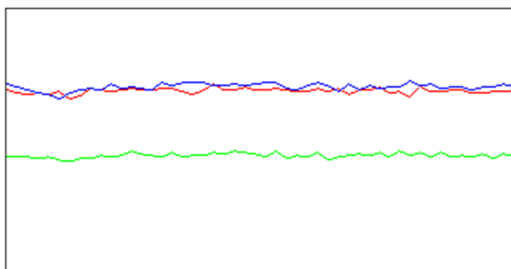

Sampling Period [sec] 0.311

Integration Time [sec] 0.1

Mass	Range	Count	Resp [cps/ug/l]	Resp (Required) [cps/ug/l]	Resp (Flag)
7	20000	13829	138288.65	34000.00	
89	50000	22287	222870.08	130000.00	
205	20000	14126	141255.44	48000.00	

Mass	Resp Ratio	Resp Ratio (Required)	Resp Ratio (Flag)
7		-	
89		-	
205		-	

Mass	RSD%	RSD% (Required)	RSD% (Flag)
7	1.483	10.000	
89	2.145	10.000	
205	1.917	10.000	

Mass	Background	Background (Required)	Background (Flag)
7	0.000		
89	0.000		
205	0.000		

Oxide/Doubly Charged Ratio


Oxide 156 / 140 0.790 %

Doubly Charged 70 / 140 2.471 %

Resolution/Axis

Q2

Tune Check Report

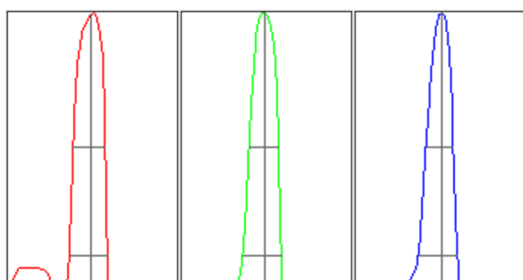

Integration Time [sec] 0.1
Acquisition Time [sec] 22.74
Y Axis Linear

Mass	Peak Height	Axis	Axis (Required)	Axis (Flag)
7	13158.06	7.00	6.90 - 7.10	
89	21374.64	89.00	88.90 - 89.10	
205	13757.27	205.05	204.90 - 205.10	

Mass	W-50%	W-10%	W-10% (Required)	W-10% (Flag)
7	0.58	0.677	0.800	
89	0.53	0.683	0.800	
205	0.49	0.691	0.800	

Tune Parameters

Plasma Parameters

Plasma Mode	---	Nebulizer Gas	0.80 L/min	Makeup Gas	0.30 L/min
RF Power	1550 W	Option Gas	0.0 %	Auxiliary Gas	0.90 L/min
RF Matching	1.30 V	Nebulizer Pump	0.10 rps	Plasma Gas	15.0 L/min
Sample Depth	8.0 mm	S/C Temp	2 °C		

Lens Parameters

Extract 1	-11.4 V	Q1 Entrance	-50.0 V	Cell Exit	-50 V
Extract 2	-200.0 V	Q1 Exit	-1.0 V	Deflect	14.6 V
Omega Bias	-95 V	Cell Focus	0.0 V	Plate Bias	-50 V
Omega Lens	8.6 V	Cell Entrance	-40 V		

Q1 Parameters

Q1 Mass Gain	131	Q1 Bias	-8.0 V
Q1 Mass Offset	137	Q1 Prefilter Bias	-6.5 V
Q1 Axis Gain	1.0011	Q1 Postfilter Bias	-3.0 V
Q1 Axis Offset	-0.21		

Q1 Ion Guide Parameters

SLS Factor	0.60	SLG Factor	0.90
------------	------	------------	------

Cell Parameters

Use Gas	No	3rd Gas Flow	0 %	Axial Acceleration	0.0 V
He Flow	0.0 mL/min	4th Gas Flow	0 %	OctP RF	170 V
H2 Flow	0.0 mL/min	OctP Bias	-8.0 V	Energy Discrimination	5.0 V

Q2 Parameters

Q2 Mass Gain	130	Q2 Axis Gain	1.0006	Q2 Bias	-3.0 V
Q2 Mass Offset	126	Q2 Axis Offset	-0.04		

Q2 Ion Guide Parameters

SLS Factor	0.50	SLG Factor	1.00
------------	------	------------	------

Wait Time Offset

Wait Time Offset	0 msec
------------------	--------

Tune Check Report

Hardware Settings

Torch


Torch H 0.4 mm Torch V -0.1 mm

EM

Discriminator 4.8 mV Analog HV 2425 V Pulse HV 1350 V

[H2]

Sensitivity

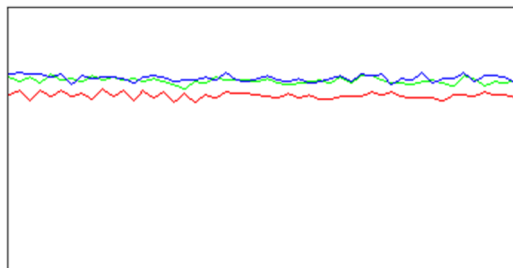

Sampling Period [sec] 0.31

Integration Time [sec] 0.1

Mass	Range	Count	Resp [cps/ug/l]	Resp (Required) [cps/ug/l]	Resp (Flag)
59 -> 59	20000	13391	133913.22		
89 -> 89	20000	14531	145312.85	90000.00	
205 -> 205	10000	7363	73633.08		

Mass	Resp Ratio	Resp Ratio (Required)	Resp Ratio (Flag)
59 -> 59		-	
89 -> 89		-	
205 -> 205		-	

Mass	RSD%	RSD% (Required)	RSD% (Flag)
59 -> 59	1.795		
89 -> 89	1.589	10.000	
205 -> 205	1.677		

Mass	Background	Background (Required)	Background (Flag)
59 -> 59	0.100		
89 -> 89	0.100		
205 -> 205	0.000		

Oxide/Doubly Charged Ratio

Oxide 156 -> 156 / 140 -> 140 0.765 %

Doubly Charged 70 -> 70 / 140 -> 140 3.049 %

Tune Parameters

Plasma Parameters

Plasma Mode	---	Nebulizer Gas	0.80 L/min	Makeup Gas	0.30 L/min
RF Power	1550 W	Option Gas	0.0 %	Auxiliary Gas	0.90 L/min
RF Matching	1.30 V	Nebulizer Pump	0.10 rps	Plasma Gas	15.0 L/min
Sample Depth	8.0 mm	S/C Temp	2 °C		

Lens Parameters

Extract 1	-9.5 V	Q1 Entrance	-50.0 V	Cell Exit	-70 V
Extract 2	-205.0 V	Q1 Exit	-3.0 V	Deflect	8.0 V
Omega Bias	-105 V	Cell Focus	4.0 V	Plate Bias	-50 V

Tune Check Report

Omega Lens 7.4 V Cell Entrance -50 V

Q1 Parameters

Q1 Mass Gain 131 Q1 Bias -3.0 V
Q1 Mass Offset 137 Q1 Prefilter Bias -8.5 V
Q1 Axis Gain 1.0011 Q1 Postfilter Bias -7.0 V
Q1 Axis Offset -0.21

Q1 Ion Guide Parameters

SLS Factor 0.60 SLG Factor 0.90

Cell Parameters

Use Gas Yes 3rd Gas Flow 0 % Axial Acceleration 2.0 V
He Flow 0.0 mL/min 4th Gas Flow 0 % OctP RF 180 V
H2 Flow 3.5 mL/min OctP Bias -1.9 V Energy Discrimination -8.0 V

Q2 Parameters

Q2 Mass Gain 130 Q2 Axis Gain 1.0006 Q2 Bias -9.9 V
Q2 Mass Offset 126 Q2 Axis Offset -0.04

Q2 Ion Guide Parameters

SLS Factor 0.50 SLG Factor 1.00

Wait Time Offset

Wait Time Offset 0 msec

Hardware Settings

Torch


Torch H 0.4 mm Torch V -0.1 mm

EM

Discriminator 4.8 mV Analog HV 2425 V Pulse HV 1350 V

[O2]

Sensitivity

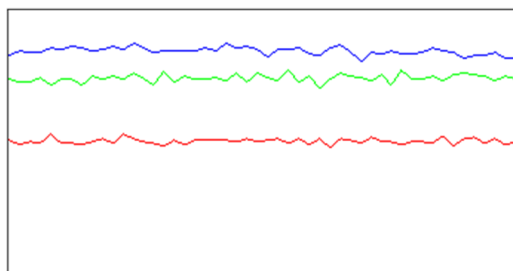

Sampling Period [sec] 0.309

Integration Time [sec] 0.1

Mass	Range	Count	Resp [cps/ug/l]	Resp (Required) [cps/ug/l]	Resp (Flag)
59 -> 59	20000	10092	100915.68		
89 -> 105	20000	14868	148684.82	70000.00	
205 -> 205	10000	8463	84634.54		
Mass	Resp Ratio	Resp Ratio (Required)	Resp Ratio (Flag)		
59 -> 59		-			
89 -> 105		-			
205 -> 205		-			

Tune Check Report

Mass	RSD%	RSD% (Required)	RSD% (Flag)
59 -> 59	2.194		
89 -> 105	1.957	10.000	
205 -> 205	1.687		

Mass	Background	Background (Required)	Background (Flag)
59 -> 59	0.000		
89 -> 105	0.000		
205 -> 205	0.000		

Tune Parameters

Plasma Parameters

Plasma Mode	---	Nebulizer Gas	0.80 L/min	Makeup Gas	0.30 L/min
RF Power	1550 W	Option Gas	0.0 %	Auxiliary Gas	0.90 L/min
RF Matching	1.30 V	Nebulizer Pump	0.10 rps	Plasma Gas	15.0 L/min
Sample Depth	8.0 mm	S/C Temp	2 °C		

Lens Parameters

Extract 1	-12.0 V	Q1 Entrance	-50.0 V	Cell Exit	-50 V
Extract 2	-210.0 V	Q1 Exit	-1.0 V	Deflect	5.4 V
Omega Bias	-100 V	Cell Focus	4.0 V	Plate Bias	-50 V
Omega Lens	8.6 V	Cell Entrance	-40 V		

Q1 Parameters

Q1 Mass Gain	131	Q1 Bias	-2.0 V
Q1 Mass Offset	137	Q1 Prefilter Bias	-6.5 V
Q1 Axis Gain	1.0011	Q1 Postfilter Bias	-9.0 V
Q1 Axis Offset	-0.21		

Q1 Ion Guide Parameters

SLS Factor	0.60	SLG Factor	0.90
------------	------	------------	------

Cell Parameters

Use Gas	Yes	3rd Gas Flow	0 %	Axial Acceleration	1.5 V
He Flow	0.0 mL/min	4th Gas Flow	25 %	OctP RF	180 V
H2 Flow	0.0 mL/min	OctP Bias	-2.5 V	Energy Discrimination	-8.0 V

Q2 Parameters

Q2 Mass Gain	130	Q2 Axis Gain	1.0006	Q2 Bias	-10.5 V
Q2 Mass Offset	126	Q2 Axis Offset	-0.04		

Q2 Ion Guide Parameters

SLS Factor	0.50	SLG Factor	1.00
------------	------	------------	------

Wait Time Offset

Wait Time Offset	2 msec
------------------	--------

Hardware Settings

Torch

Torch H	0.4 mm	Torch V	-0.1 mm
---------	--------	---------	---------

EM

Discriminator	4.8 mV	Analog HV	2425 V	Pulse HV	1350 V
---------------	--------	-----------	--------	----------	--------